



**Instructions For Use
Medial Relief Brim Installation**

©2017 LIM Innovations MKT1.0-111617

MEDIAL RELIEF IS ON THE WAY.

YOU SPOKE. WE LISTENED.

This document has been developed to support the implementation of the medial relief brim to the Infinite Socket™ TF system. It is intended for use by clinicians and industry professionals, and will provide benefits in establishing an economic and effective approach to the control of clinical outcomes.

MORE QUESTIONS?

CONTACT OUR
CUSTOMER SERVICE
(844) 888-8546




424 9th Street
San Francisco, CA
94103

844.888.8546
LIMinnovations.com



NEW COMPONENTS

- A. Ischial Seat Base
- B. Threaded Washer
- C. M5x16mm Socket Head Screw*
- D. Medial Relief Brim (MRB)
- E. TF Ischial Strut with integrated inserts

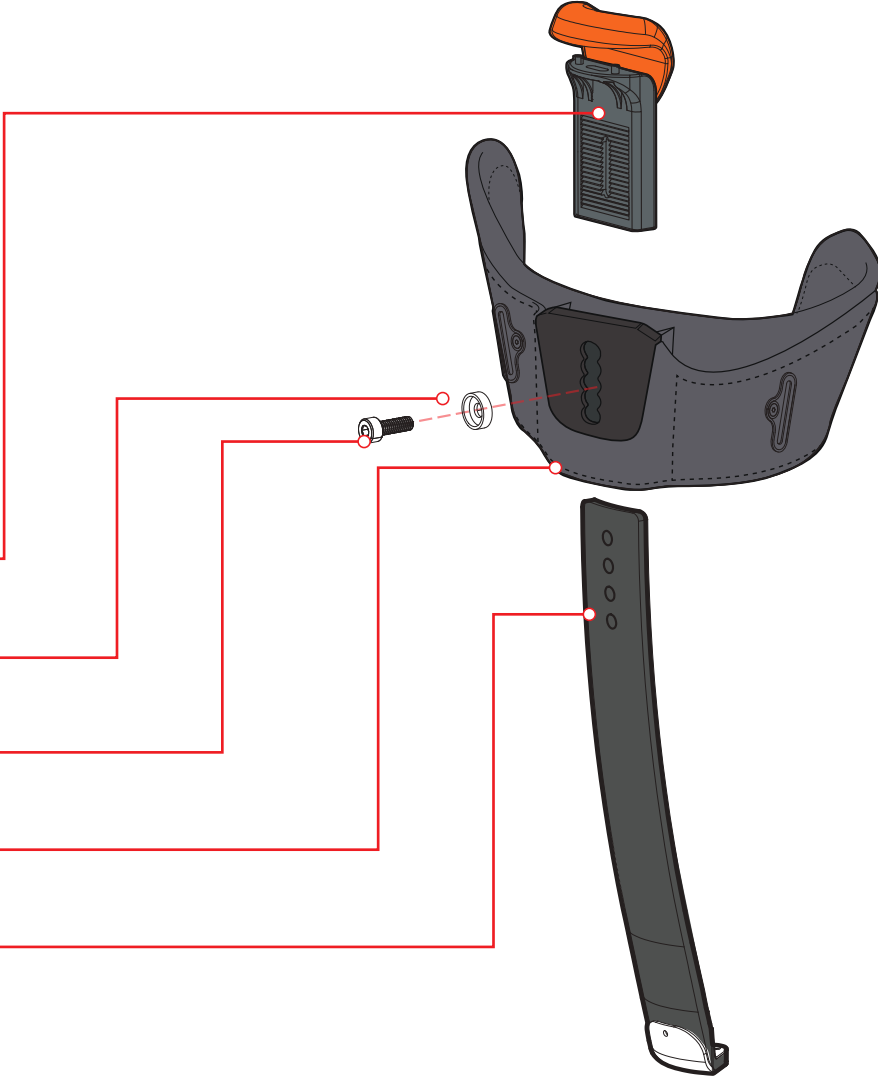

*Ischial seat screw torque setting: 3 N.m



FEATURES

The medial relief brim allows for 3 adjustments.

Ischial seat and medial relief brim macro adjustment
4 holes on the strut, 15 mm apart

Ischial seat height micro adjustment
-Micro teeth allow for an adjustment of 30 mm total

Medial relief brim height micro adjustment
-5 holes on the brim base, 8.8 mm apart

